



- ▶ Automated needle retraction
- ▶ Effectively reduces the risk of needlestick injuries
- ▶ Quick and easy to use
- ▶ Activation can be achieved in a variety of ways
- ▶ May be activated before or after removal from patient
- ▶ Pre-removal activation reduces the risk of exposure to the contaminated sharp
- ▶ Activation of safety mechanism automatically covers syringe luer tip, reducing risk of exposure to medication, blood/body fluids
- ▶ Remains safe through disposal





compatible with
both luer lock &
luer slip syringes

may be used with
pre-filled syringes



(50 per box /400 per case)

NEEDLES ONLY			
Catalog #	Needle Size	Catalog #	Needle Size
For syringes 3mL and smaller		For syringes 10mL and smaller	
82211-USA	30G x 1/2"	85291-USA	25G x 5/8"
82091-USA	25G x 5/8"	85211-USA	25G x 1"
82011-USA	25G x 1"	85011-USA	25G x 1 1/2"
82031-USA	23G x 1"	85231-USA	23G x 1"
		85021-USA	22G x 1"
		85031-USA	22G x 1 1/2"
		85041-USA	21G x 1"
		85051-USA	21G x 1 1/2"
		85061-USA	20G x 1"
		85071-USA	20G x 1 1/2"
		85081-USA	18G x 1"
		85091-USA	18G x 1 1/2"

WITH 3mL SYRINGES		
	Catalog #	Needle Size
	82391-USA	25G x 5/8"
	82311-USA	25G x 1"
	85311-USA	25G x 1 1/2"
	82331-USA	23G x 1"
	85321-USA	22G x 1"
	85331-USA	22G x 1 1/2"
	85341-USA	21G x 1"
	85351-USA	21G x 1 1/2"
	85361-USA	20G x 1"
	85371-USA	20G x 1 1/2"
	85381-USA	18G x 1"
	85391-USA	18G x 1 1/2"

OTHER AVAILABLE SAFETY PRODUCTS

- ▶ VanishPoint® Syringe
- ▶ VanishPoint® Blood Collection Tube Holder
- ▶ VanishPoint® Blood Collection Set
- ▶ VanishPoint® IV Catheter
- ▶ EasyPoint® Blood Collection Tube Holder with Needle



RETRACTABLE
TECHNOLOGIES, INC.®

To order these products,
please contact

rtisales@retractable.com

511 Lobo Lane • Little Elm
Texas • 75068-5295 • USA
Tel.: (972) 294-1770
Fax: (972) 294-4400
Toll Free: 1-888-703-1010

© 2015, 2025 Retractable Technologies, Inc. All rights reserved. MKT-BRO-EPN-0002 REV.0 EN, AUG 2025



SAVES MONEY, TIME, AND LIVES



www.retractable.com




RETRACTABLE
TECHNOLOGIES, INC.®

easypoint™

RETRACTABLE NEEDLE

SAVES MONEY, TIME, AND LIVES

AUTOMATED
RETRACTION
- now with -
ADDED
VERSATILITY

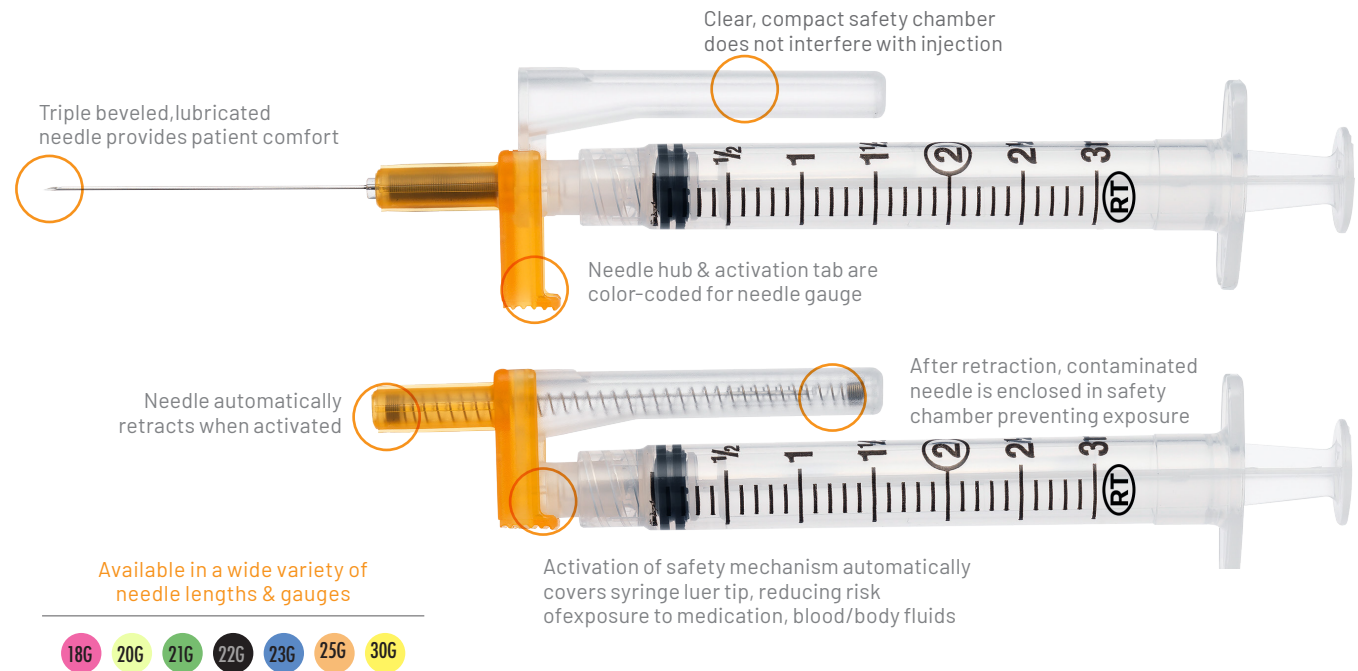

easy to use


- Minimal change in injection technique
- Allows for single-handed activation
- Clinician's hands and fingers remain behind the exposed sharp at all times

1



Securely attach EasyPoint needle to a syringe. Administer medication according to institutional policy.

2

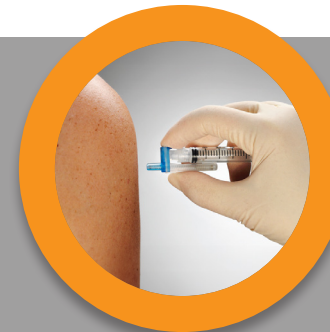

Apply pressure to the color-coded tab and the safety chamber using a squeezing motion, allowing the needle to automatically retract into the safety chamber.

3



Dispose in an appropriate sharps container.