




- ▶ Automated needle retraction
- ▶ Effectively reduces the risk of needlestick injuries
- ▶ Quick and easy to use
- ▶ Activation can be achieved in a variety of ways
- ▶ May be activated before or after removal from patient
- ▶ Pre-removal activation reduces the risk of exposure to the contaminated sharp
- ▶ Activation of safety mechanism automatically covers syringe luer tip, reducing risk of exposure to medication, blood/body fluids
- ▶ Remains safe through disposal

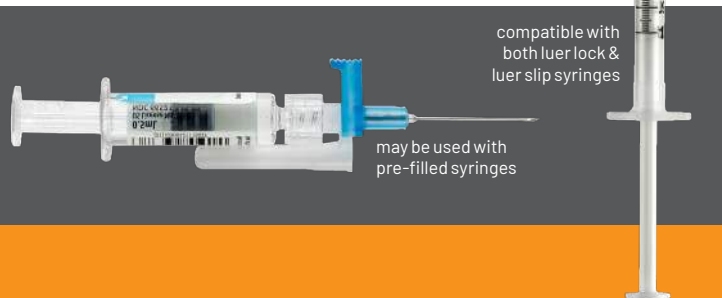

(50 per box / 400 per case)

NEEDLES ONLY						
Catalog #	Needle Size	Catalog #	Needle Size			
For syringes 3mL and smaller		For syringes 10mL and smaller				
82211	30G x 1/2"	85291	25G x 5/8"		82391	25G x 5/8"
82091	25G x 5/8"	85211	25G x 1"		82311	25G x 1"
82011	25G x 1"	85011	25G x 1 1/2"		85311	25G x 1 1/2"
82031	23G x 1"	85231	23G x 1"		82331	23G x 1"
		85021	22G x 1"		85321	22G x 1"
		85031	22G x 1 1/2"		85331	22G x 1 1/2"
		85041	21G x 1"		85341	21G x 1"
		85051	21G x 1 1/2"		85351	21G x 1 1/2"
		85061	20G x 1"		85361	20G x 1"
		85071	20G x 1 1/2"		85371	20G x 1 1/2"
		85081	18G x 1"		85381	18G x 1"
		85091	18G x 1 1/2"		85391	18G x 1 1/2"

OTHER AVAILABLE SAFETY PRODUCTS

- ▶ VanishPoint® Syringe
- ▶ VanishPoint® Blood Collection Tube Holder
- ▶ VanishPoint® Blood Collection Set
- ▶ VanishPoint® IV Catheter
- ▶ EasyPoint® Blood Collection Tube Holder with Needle



RETRACTABLE
TECHNOLOGIES, INC.®

To order these products,
please contact

rtisales@retractable.com

511 Lobo Lane • Little Elm
Texas • 75068-5295 • USA
Tel.: (972) 294-1010
Fax: (972) 294-4400
Toll Free: 1-888-703-1010



SAVES MONEY, TIME, AND LIVES



www.retractable.com



RETRACTABLE
TECHNOLOGIES, INC.®



SAVES MONEY, TIME, AND LIVES

**AUTOMATED
RETRACTION**
- now with -
**ADDED
VERSATILITY**

Triple beveled,
lubricated needle
provides patient
comfort

Needle hub & activation
tab are color-coded for
needle gauge

Clear, compact safety chamber
does not interfere with injection

Needle automatically
retracts when activated

Activation of safety mechanism automatically
covers syringe luer tip, reducing risk of
exposure to medication, blood/body fluids

After retraction, contaminated
needle is enclosed in safety
chamber preventing exposure

Available in a wide variety of
needle lengths & gauges



easy to use

- Minimal change in injection technique
- Allows for single-handed activation
- Clinician's hands and fingers remain behind the exposed sharp at all times



Securely attach EasyPoint
needle to a syringe. Administer
medication according to
institutional policy.



Apply pressure to the color-coded
tab and the safety chamber using
a squeezing motion, allowing the
needle to automatically retract into
the safety chamber.



Dispose in an appropriate
sharps container.